

READY-TO-USE PRODUCTS

SURFA'SAFE PREMIUM



- » Disinfectant detergent formulated without alcohol, without perfume nor CMR.
- » Quick antimicrobial efficacy: bacteria including BHR, Rotavirus... from 3 min.
- » Quick drying, does not leave any trace.
- » Wide compatibility with materials, specially polymers (polycarbonates...).



INDICATIONS

Detergent disinfectant foam, for cleaning and disinfection of surfaces, non-immersible and non-invasive medical devices (conventional ultrasound probes, cables and connectors, pressure sensors, blood sugar testers, tensiometers, reflex hammers...) and external structures of equipments (respirators, haemodialysis generators...).

COMPOSITION

Didecyltrimethylammonium chloride, excipients.

PACKAGING - REF.


12 bottles, 750ml each, with foaming dispensers.....2419.544
 16 bottles, 500ml each2419.105

INSTRUCTIONS FOR USE



READY-TO-USE
Active from 3 min.

ACTIVE AGAINST	ACCORDING TO STANDARDS
BACTERIA	EN 1040, EN 13727, EN 14561, EN 13697 Active against MDRB and HRB according to EN 13727, EN 13697
MYCOBACTERIA	Mycobacterium terrae : EN 14348, EN 13697
YEASTS/MOULDS	EN 1275, EN 13624, EN 13697, EN 14562
VIRUSES	Adénovirus and Norovirus according to EN 14476 According to EN 14476 : PRV (surrogate of HBV), BVDV (surrogate of HCV), Vaccinia virus, Herpesvirus, Rotavirus, Polyomavirus SV40

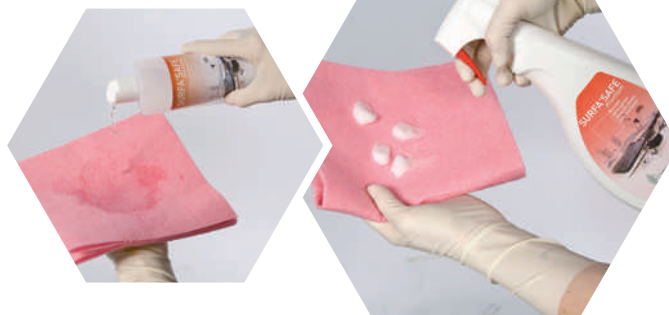

Exists in wipes

WIP'ANIOS EXCEL



- » SURFA'SAFE PREMIUM impregnation solution.
- » Formulated without CMR substances, without perfume, nor alcohol.
- » 100% biodegradable viscose wipe.



INDICATIONS

WIP'ANIOS EXCEL wipes are impregnated with a cleaning and disinfecting solution for the treatment of surfaces and medical devices.

- **Cleaning and disinfection of non-immersible and non invasive medical devices** (ultrasound probes, cables and connectors, pressure sensors, blood sugar testers, tensiometers...) and external structures of equipments (respirators, haemodialysis generators...).
- **Wiping of endoscopy sheaths** during preliminary treatment, prior to immersion for cleaning and disinfection.

PACKAGING - REF.

6 sachets, 100 wipes each..... 2446.655

INSTRUCTIONS FOR USE



SINGLE USE WIPES

Active from 3 min.



MORE DETAILS IN THE BROCHURE AND VIDEO

Our advice

- For large surfaces, preferably use the 500ml bottle.